

LEAD-LINED

UNIT DOSE CABINET

The Unit Dose Cabinet is designed for hot labs with limited space.

The upper left compartment holds two Sharps containers to facilitate decay rotation. A small chute with removable shielded lid allows syringe disposal into the front container. A large rear port allows safe removal of the second decayed container. The upper right compartment is for bulk storage of unit dose ammo boxes. The shielded door provides access without additional exposure.

Side-by-side middle drawers can be used for storage of syringes, sources, pigs, radioisotopes and other small items requiring lead shielding.

Phantoms and flood sources can be stored in the bottom drawer which runs the full width of the cabinet.

SPECIFICATIONS:

Dimensions: 30.5" w x 24" depth x 36.5" h (77.5 x 61 x 92.7 cm)

Lead Shielding: 25" (.64 cm), .5" (1.3 cm) or 1" (2.5 cm) thick on all six sides

Drawer Dimensions:

Bottom Drawer: I.D.: 23.5" w x 21.8" depth x 4.5" h (59.7 x 55.4 x 11.4 cm)

Left & Right Drawers: I.D.: 9.38" w x 17.88" depth x 4.5" h (23.8 x 45.4 x 11.4 cm)

Upper Doors: I.D.: 11.5" w x 17.62" depth x 11" h (29.2 x 44.7 x 27.9 cm)

Doors and Drawers: Key-locked

Countertop: Stainless steel with 4" backplash and .5" spill-proof lip


Finish: Durable powder-coat paint

Weight:

244-120: 786 lb (357 kg)

244-121: 1293 lb (586 kg)


244-122: 2328 lb (1056 kg)

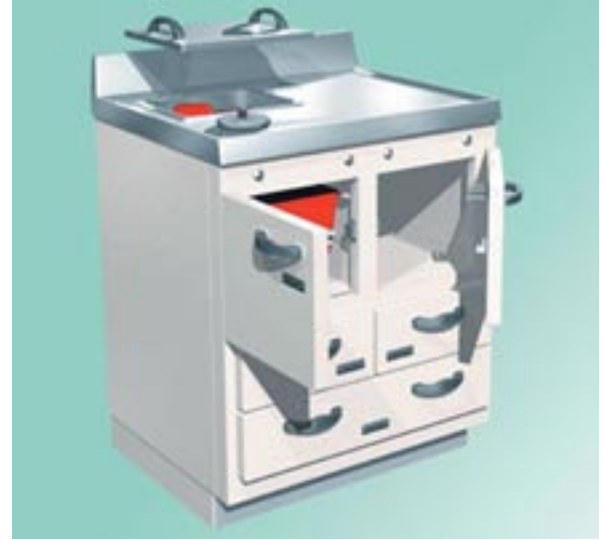

244-120 Cabinet, Unit Dose, .25" lead ... \$12,100.00

244-121 Cabinet, Unit Dose, .5" lead 13,600.00

244-122 Cabinet, Unit Dose, 1" lead 16,400.00

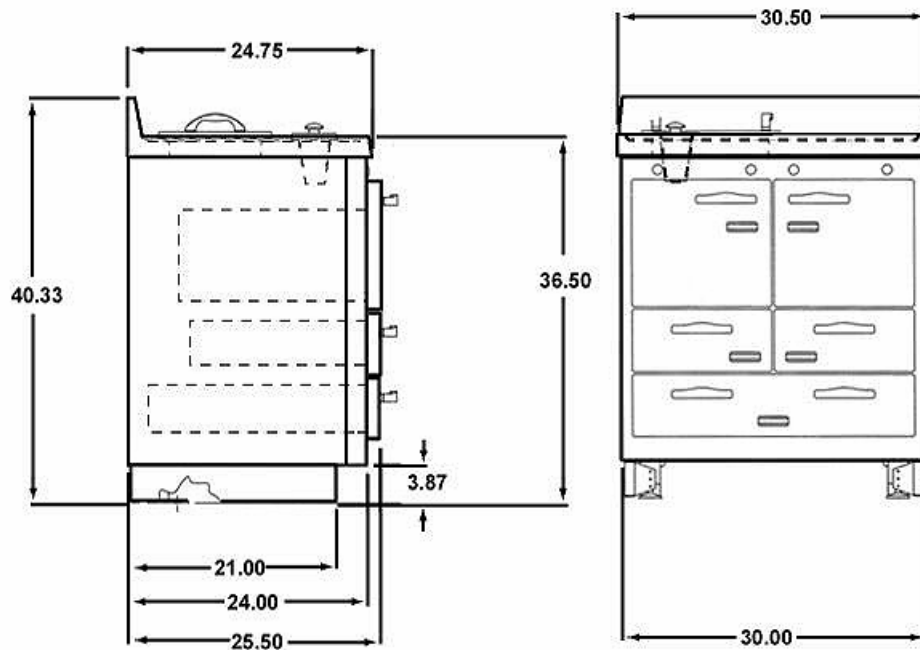

Prices subject to change without notice.

www.jgravengard.com 2007-10-12

J. Gravengard Corporation

Medical Imaging Sales & Service

PO Box 20171 Portland, Oregon 97294 Nationwide (800)796-SCAN

Local (503)255-0535 Fax (503)255-7047